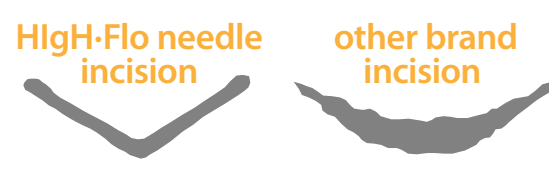
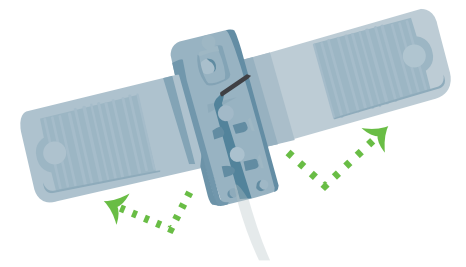
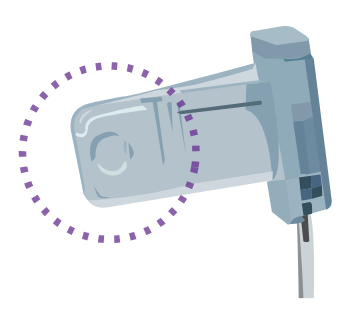
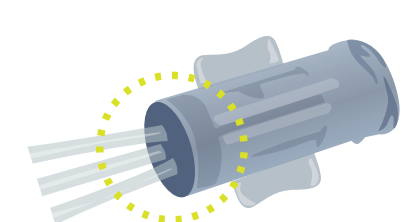
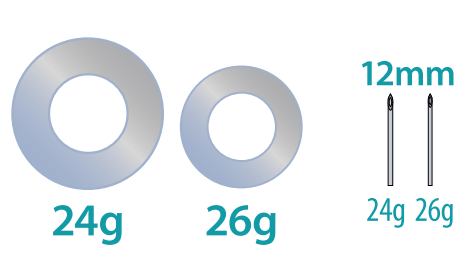
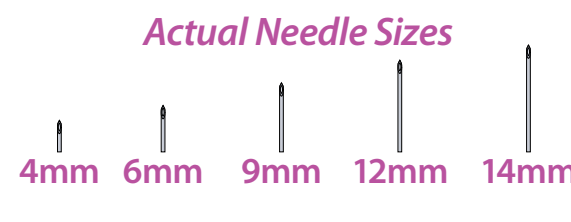
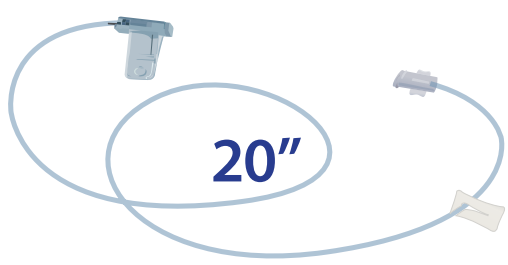

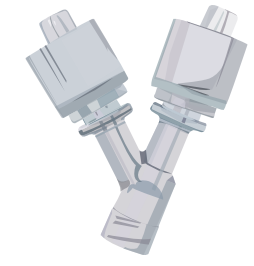



get to know


High Flo™

Subcutaneous Safety Needle Sets

What makes a High Flo needle so unique?

Needle Tip Design	5-bevel backcut needle tip design minimizes damage, scarring and local site reactions.	
Flying Hinge	Flying hinge butterfly wings have great flexibility and allow each needle to stay in place.	
Safety Closure	Safety closure protects against accidental needle sticks.	
Custom Luer	Custom luer design provides even flow to every site.	
Gauge Options	26g needles designed for sub-q flow rates up to 40ml/hr/site and 24g for up to 300ml/hr/site.	
Length Options	Needles are available in 4mm, 6mm, 9mm, 12mm, and 14mm lengths.	
Tubing Length	Efficient and practical 20" length (in 26g sets) for less tangles and faster flow rates. 24" in 24g sets.	
Extension Set	Extension set available for those who prefer a longer length.	
Y-Connector	Combines sets for more than 6 sites (up to 8), or to enable infusing with different length needles.	
3M Tegaderm™	Genuine 3M Tegaderm™ transparent dressing is included with all sets.	

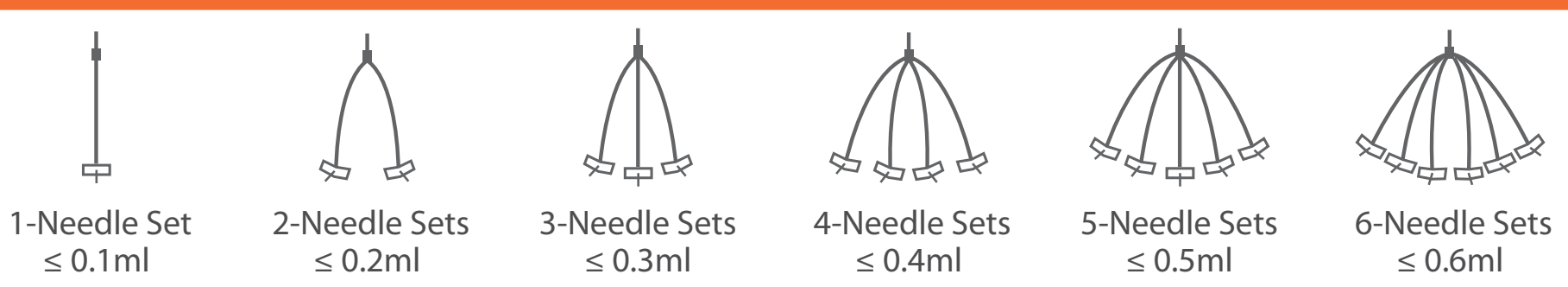26g High Flo Needle Sets Options (Residual Volume)



24g High Flo Needle Sets Options (Residual Volume)

