

CFUR30I & CFUR40I


Urology Tables

- 3 table extensions included
- 500 lb Patient Weight Capacity
- Warranty: 3 years parts, 2 years labor



FM584901

made in the USA
with US & imported parts

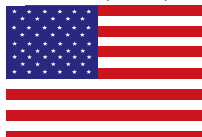

OAKWORKS[®]
MEDICAL PRODUCTS

3 Radiolucent Extensions

33" WIDE SUPPORT SURFACE

The two lateral extensions give you the support area you need for patients' arms while maximizing access to the 23" W x 47" L Carbon Fiber Top. Imaging Area 19" W x 24" L (48cm x 61cm)

Radiolucent patient transfer extension creates 51" long metal-free imaging area between the Trails.



(2) 5" X 24" RADIOLUCENT LATERAL EXTENSIONS #81974



Removable lateral extensions slide and lock on T-rails.

(1) 23" WIDE X 28" LONG RADIOLUCENT PATIENT TRANSFER EXTENSION #81972



Lightweight design can be used for general imaging when a full length diving board table design is required.

CFUR301 #81818

CFUR401 #81817

STANDARD	
Carbon Fiber Top 53" (135cm) Metal free imaging space with Patient Transfer Extension	23" (58cm) wide x 47" (119cm) long
CFUR301 MOVEMENTS	
Motion 1	Electric Height Range 25.5" - 43.5" (65-110cm)
Motion 2	Electric Trendelenburg tilt +15°/-11°
Motion 3	Electric Longitudinal travel 10" (25cm)
CFUR401 MOVEMENTS	
Motion 1	Electric Height Range 26.5" - 44.5" (67-114cm)
Motion 2	Electric Trendelenburg tilt +15°/-11°
Motion 3	Electric Longitudinal travel 10" (25cm)
Motion 4	Electric Lateral travel ±4" (10cm)
Controls	
Controls	Multi-function hand control
Weight Capacity	
Weight Capacity	500 lb (227kg) Each Extension: 200 lb (91kg)
Product Weight	
Product Weight	CFUR301 568 lb (258kg) CFUR401 590 lb (268kg)
Radiolucent Lateral Extension	
Radiolucent Lateral Extension	(2) 5" x 24" (12cm x 61cm)
Radiolucent Patient Transfer Extension	
Radiolucent Patient Transfer Extension	(1) 23" x 28" (58cm x 71cm)
Included Features	
Included Features	1" (3cm) Radiolucent Comfort Foam™ table top pad, 4" (10cm) locking casters, Digital angle gauge, patient restraint straps (2), Industry standard full length T-rails on each side, Urology bag hoop
Warranty	
Warranty	3 years parts, 2 years labor
Safety Listings	
Safety Listings	FDA listed, CE marked, cETLus listed, CB Certificate

OPTIONS must be ordered with table.		#
5" (13cm) casters	upgrade from the standard 4" (10cm)	65427
Foot Control (CFUR301)	multi-functional	75076
Foot Control (CFUR401)	multi-functional	75078
Non-USA electric system	please contact us for order details	73134
ACCESSORIES order any time.		#
Protective Drape	Clear PVC - attaches to T-Rail	81961
Radiolucent Lateral Extension	5" x 24" (13cm x 61cm)	81962
Radiolucent Lateral Extension	5" x 14" (13cm x 36cm)	69794

UROLOGY BAG HOOP

Table comes with hoop and clamps to t-rails. (bags not included)

INCLUDED



5" CASTER #65427

Recommended for mobile service providers

OPTIONAL





Protective Drape #69390

INCLUDED FEATURES:

(3) Radiolucent Extensions

(2) Patient Restraint Straps

84" L with 24" of Velcro® at each end



1" Comfort Foam™ Table Top Pad



4" Locking Casters


Industry Standard Full Length T-rails on each side

The Integrated T-rails are US standard size of .375" (9.5mm) thick by 1.25" (28.6mm) wide.

Digital Angle Gauge

Urology Bag Hoop

Model: Urocatch® with intergrak clamps
(Bag not included)

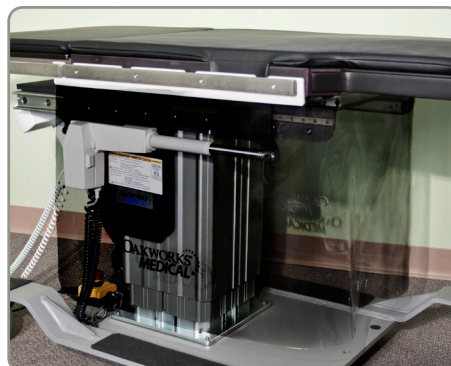

This removable clear vinyl drape acts as a splash shield. It attaches to the t-rail and facilitates movement of C arms and lithotripters without concerns about bumping into the table mechanisms.



FM584901

Oakworks Medical is ISO 13485:2003 certified. The ISO standard is intended to ensure that medical device manufacturers have the required quality management systems in place to safely design, manufacture and distribute medical devices. This registration certifies that Oakworks Medical's quality system conforms to these requirements in our state-of-the-art 90,000 square foot factory in Pennsylvania.

717.759.3146 www.oakworksmed.com email: customerservice@oakworks.com